## Vibramat

Massagerät | portable massage

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**EN** Operation manual

Deutsch.....	3
English.....	15

# Deutsch

## Inhalt

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# English

## Table of Contents

Before use .....	16
Safety precautions .....	16
Technical data .....	16
Basic equipment .....	17
Optional .....	17
Standard equipment .....	17
Special accessories .....	17
Detailed massage instructions .....	18
Optimal effects of a complete massage .....	18
The massage should be in counter-clockwise direction .....	19
Up to 10 massage variations .....	19
Start with your right foot .....	19
Arm and leg massage .....	20
Unwanted inches .....	21
Connective tissue massage .....	21
Neck, Back .....	21
Muscle massage .....	22
Sport massage .....	22
Local massage .....	22
Internal organs .....	23
Bronchial asthma and emphysema .....	23
Lumbago .....	23
Special care .....	23
Accessories' upkeep: .....	23
Chart Massage Schedule .....	24
Guarantee .....	25
General Terms and Conditions of Service .....	26

## Before use

No massage in case of feverish colds, varicose veins, thrombosis, phlebitis, pregnancy, jaundice, tumours or other acute ignitions. In case of any kind of illness you always have to consult your physician before starting the massage, when practising special application forms always await the prescription of your physician.

The massage appliance is a product which can help you to enhance and promote your well-being. It is not a medical product according to the medicine-product-law.

Please read this instruction booklet carefully and keep it handy for future reference.

## Safety precautions

This appliance conforms to the recognized technical specifications on safety and corresponding safety regulations.

It must only be connected to an AC power supply of the same voltage as that shown in unit identification plate.

Never pull the plug out of its socket by the lead.

Always pull out the mains plug before cleaning or working on the appliance.

This equipment must be disconnected from the mains when not in use.

Children must be supervised during use.

Please do not cover the air intake while using. Overheating-Danger!

This device has not been designed for being used by persons (including children) that

are physically or mentally handicapped or have sensory disabilities. Moreover, it is not allowed to use this device without sufficient experience and/or knowledge, unless these persons will be supervised by persons responsible for their security or in case they have been instructed how to use this device.

Children are to be supervised in order to make sure that they do not play with this device.

### ATTENTION!

Do not plunge the massager into water or expose it to high humidity.  
Risk of electric shock! Danger to life!

### CAUTION!

Do not plunge the massager into water or expose it to high humidity.  
Risk of electric shock! Danger to life!

## Specifications

Power supply: 230V 1N AC 50Hz

Power output: 100 W

Motor: special high-performance motor, interference suppressed, completely insulated, maintenance free. Suitable for long-term use.

Speed: 1400 - 2800 rpm

Operation: 2-speed switch for mild and strong massage.

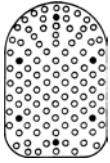
Casing and base plate made of shock-resistant plastic.

Mains lead: CEE-Norm plug, 2.8 m mains cable.

This unit complies with the 82/499/EC requirements.


## Basic equipment

### Base plate

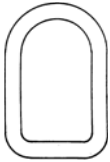
pre-mouled

### Cellular rubber pad

Item. 2000 7364

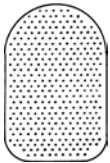Included in delivery


### Plastic hyperemization brush

Item 2000.7365

Included in delivery
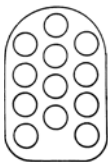
## Optional accessories

### Studded cellular rubber shoe

Item. 2000.7226

## Standard features

### Massage pads and brushes

The **plastic hyperemization brush** and the **cellular rubber pad** are included with the massager. The **studded cellular rubber shoe** can be purchased separately.

If you wish to obtain better blood circulation use, according to your sensitivity, one of the two hyperemia brushes.

Hyperamia means: The process by which you improve the blood circulation. The cellular rubber pad is more gentle on your skin and is softer than the plastic hyperemization brush pad.

The smooth cellular rubber pad cushions the massage strokes of the base plate, but at the same time, increases the in-depth effect of the massage.

To put on the cellular rubber pad or the studded cellular rubber shoe, at first push the straight rear end of the base plate into the corresponding part of the shoe/pad. Then pull the rounded front part of the pad over the rounded front of the bottom plate.

The procedure to put on the nylon brush is the same. First it is important that you fit the straight rear end of the brush over the straight rear end of the base plate.

Then press the rounded part of the brush over the rounded front section of the base plate. In order to take off these shoes, press your thumbs firmly on both sides of the back of the pad and pull the pad off the base plate.

Notice: It is normal that you need to apply certain force in order to pull the shoe on the base plate. The shoes are slightly smaller than the plate to ensure a tight fit. So that the shoe cannot get off during an intensive massage.

## Massage instructions

Take your time to get yourself familiar with your new massager. The more you know about the possibilities of this device, the more you are going to benefit from it.

The high massage efficiency of this device is based on the unique their oscillatory action which combines the horizontal vibration with the vertical impulses.

This provides a 3D effect - massage action on the length, the width and into the depth. Muscles and tissue are loosened, the blood circulation is improved and the metabolism is enhanced.

Contrary to a time-consuming and tiring hand massage, the large-surface massager does the work 3 to 4 times faster and effectively covers large surfaces. It even offers another advantage - from the first application, you lose the fear of the often painful hand massage. Large-surface massagers rarely cause pain.

**No massage in the kidney area, as this may lead to kidney damages!**

The functional and ergonomic design of this device combines the proven technology with the high convenience in use and offers a maximum efficiency with continuously high performance as well as a low noise operation.

## Optimal effects of a complete massage

- **Shows: free choice\***
- **Duration of massage: 10 -12 minutes**

Do you wish to improve your general health condition?

Then a body or complete massage is suggested once a week.

While a hand massage takes 40 to 50 minutes of your free time, large-surface massager does the work in 10 to 12 minutes.

**Do not massage longer!  
A lengthy application is useless since the effects obtained could be neutralized by a longer massage.**

You can divide the complete massage on two weekdays.

On one day, for example, massage your feet, legs, posterior and arms. Three or four days later the body: chest, neck and both side of the spinal cord

\* Free means:

Either without a show at all (just the base plate) or with a hyperemization brush, a cellular rubber shoe, a studded cellular rubber show.

## The massage should be in counter-clockwise direction

Hold the massager by the handle with one hand and with the other hold the knob. Massage slowly and in circles, always in counter-clockwise direction. When massaging the whole body, follow this pattern:

1. Feet
2. Legs
3. Posterior
4. Arms
5. Chest
6. Abdomen
7. Neck
8. Back

If somebody can massage you, lie down and relax your muscles.

The massages should be given with the base plate. Press firm but gently the surface you wish to massage

**Strong pressure does not produce a stronger effect, on the contrary, it diminishes the result.**

According to your needs and the sensitivity of your skin, you may choose from the four accessories.

## Up to 10 massage variations

Do you prefer a stronger or a more mild massage effect? You may choose between the two vibration speeds.

- Speed 1 - soft
- Speed 2 - strong

And what about the massage variations?

These can be obtained by using the accessories. You can also use the massager without any accessory, just the base plate.

Furthermore you can achieve a more pointed or a surface-oriented massage by using either the flat surface or the edge of the base plate.

## Start with your right foot

- **Accessory: base plate, nylon brush or smooth rubber pad.**
- **Massage: 1- 2 minutes each foot.**

Start the body massage with your right foot.

At first massage between the two ankles and then from the heels up to the toes. Now you can massage your soles.

Massage your feet only while sitting down.

### **NOTE: Do not step on the massager!**

An additional massage from the Achilles'tendon to the calf is highly recommended.

Can you feel the results? You may now start massaging your left foot in the same way.

- Do you have to stand while working?
- Do you have swollen feet and legs?

Regular massage helps in many cases. After a few massage strokes you will realize how your tired feet become lighter. The effects of tiredness are vanishing like magic.

It is important that feet are taken care of every day since in most cases the swollen parts slowly disappear.

Flat feet can also be avoided if massaged beforehand. In case of flat feet, the massage will decrease pain, especially combined with active foot gymnastics.

The foot massage with a large-surface massager is gymnastics and massage at the same time.

### Arm and leg massage

- **Accessories: base plate, nylon brush pad.**

- **Duration: 2 -3 minutes each leg.**

The leg massage improves circulation of the blood and loosens muscular cramps in legs and thighs. Leg cramps can be treated satisfactorily. The calves, the back of your knee and the lower thighs should be also massaged. Even if you have soft skin, this type of massage will be pleasant.

**Be careful if you have varicose veins, thrombosis or phlebitis. Ask your doctor beforehand.**

Do your arms itch?

Do your hands feel numb?

Daily massage can help. Start again with your right hand. If you have heart ailment, massage your left arm vary carefully. Check with doctor first.

- **Accessories: base plate.**

- **Duration: 2 -3 minutes each arm.**

Inflexion side

Extensor side

Extensor side

Inflexion side

## Unwanted inches

- **Accessories:** see chart massage schedule.
- **Duration:** see chart massage schedule.

Place the massager where you have excessive inches - waist, thighs, hips, posterior, and stomach.

**Daily** massage will activate the tissue fluid, the retention of water will be reduced, and hence, the loss of "fat" is helped.

Massage as shown below.

For faster results go on a diet simultaneously.

## Connective tissue massage

- **Accessories:** optional.
- **Duration:** 2 - 5 minutes.

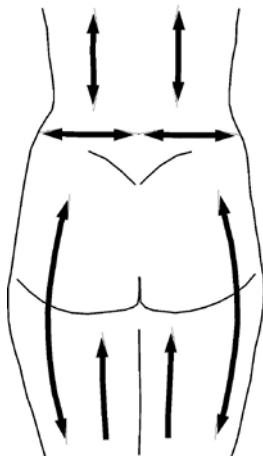
If you wish to slow down the ageing process you should start on time with the necessary measures to regenerate your cells.

A very positive way to achieve this is by a connective tissue massage.

This will bring back the tissue's elasticity.

Metabolic residues slowly disappear as well as the functional problems of the organs involved.

Massage as for the complete body massage.

## Neck, Back

- **Accessories:** see chart massage schedule.
- **Duration:** 4 minutes approximately.

Long car rides, little movement at work, too much TV watching. We sit down too much. And the results?

Strained neck and shoulder muscles. Back and neck pain.

The large surface massager is excellent for preventive spinal cord massage.

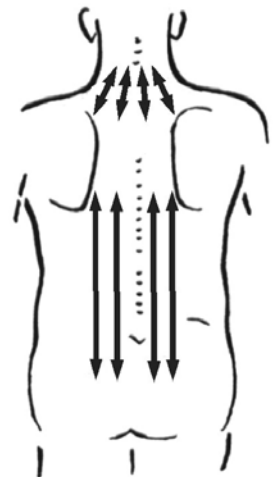
Its oscillatory movement is carried over directly to the spinal cord and cramps in the vertical back muscles will disappear.

Pains here originated can be decreased or disappear.

Massage as shown in the picture.

First the neck (approx. 1 minute). Then right and left of the spinal cord (approx. 1 minute).

**i** Damaged discs are extremely painful. Preventive vibratory massage with a large-surface massager could avoid this. Check with doctor first.

## Muscle massage

- **Accessories:** base plate, cellular rubber pad, studded rubber pad.
- **Duration:** as long as desired.

Cramped or tense muscles?

Due to its three-dimensional massage effect (oscillatory movement) and the large surface of the device you obtain optimum results.

Vibratory massage, especially on the hips, shoulders, and lower back region is very effective to loosen muscular cramps.

## Sport massage

- **Accessories:** optional.
- **Duration:** depending upon which part of the body, 1- 4 minutes.

Not only the muscles directly strained by sport are massages, but also those stimulated indirectly.

The large-surface massager is also recommended for a short message during a championship of training since not only is the massage intensive, but most important of all, it can be performed during a short period of time.

A decisive advantage:

You may massage without removing your clothes.

## Local massage

- **Accessory:** smooth rubber pad.
- **Duration:** 1- 2 minutes.

Vibratory massage not only relaxes and keeps you fit but also can be employed satisfactory for local massage.

The period of malaise can very often be shortened. Its very differential oscillating movement obtains even better results than a massage by hand.

### Warning:

Use carefully. Consult a physician before you use it for local massage.

Further information can also be obtained at book shops.

## Internal organs

- **Accessory:** cellular rubber pad.
- **Duration:** 1- 2 minutes.

Organs can be positively influenced and the metabolism improved by massaging a specific part of your body.

Ask your doctor, masseur or gym therapist for these special cases.

To stimulate bowel movement, massage your back and stomach as shown on the two pictures below. When massaging your stomach, use the slow speed.

No massage in the kidney area, as this may lead to kidney damages!

Front

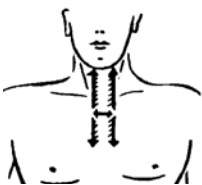
Back

## Bronchial asthma and emphysema

- **Accessories:** cellular rubber pad.
- **Duration:** according to doctor's advice.

Relief is obtained by massaging slowly and smoothly. Massage both sides of the larynx, the breastbone, and also between your spinal cord and shoulders. Check with your doctor first.

Front

Back

## Lumbago

- **Accessories:** base plate, studded rubber pad.
- **Duration:** 3 - 5 minutes.

This sudden back and neck pain needs a daily massage.

Massage the lower back at both sides of the spinal cord.

Check with your doctor first.

## Maintenance and care

This large-surface massager is absolutely maintenance free. There are no inner parts which require any servicing.

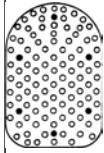
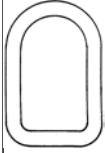
The supplementary massage shoes should be cleaned from time to time by washing them in lukewarm soapy water.

The base plate which is fastened to the massager can be cleaned with a damp cloth. DO NOT immerse the massager in water under any circumstances.

The massage shoes are wear items. Replace the broken or worn-out shoes.

## Chart Massage Schedule

Duration of single-massage, vibration speed and instructions for use

Zones:	Duration in minutes	Speed				
arms	2 - 3	I II	•			•
stomach	1 - 2	I			•	
legs	2 - 3	II	•	•		•
connective tissue	2 - 5	I II	•	•	•	•
chest	1 - 2	I			•	
double chin	up 1	I			•	
between calves and feet	1 - 2	II			•	
feet	1 - 2	II	•	•	•	
body*	10 - 12	I II	•	•	•	•
posterior	2 - 3	II	•	•		•
face	up 1	I			•	
hands	1 - 2	I II	•			•
lumbago	3 - 5	I II	•			
hips	2 - 3	II	•	•	•	
organs	1 - 2	I			•	
lumbar region	1 - 2	I			•	
muscles	1 - 5	I II	•		•	
neck	1	I II	•		•	
thighs	2 - 3	II	•	•	•	•
back	2	II				•
shoulders	2	II	•			
waist	1 - 2	II		•	•	•
calves	2 - 3	II	•	•		
spinal cord relax	1	I II			•	

\*does not equal individual massage duration for each part of the body a full massage can be performed quicker.

## Recycling

Devices or lighting elements that will not be used any longer have to be handed in at a recycling station according to regulation 2012/19/EU. Do not dispose it with the normal household waste.

Please retain this address together with the installation guide for further references.

To help us answer your questions quickly and competently please provide the information printed on the type shield including the model, item no. and serial no., in all inquiries.

**Equipment commissioning date:**

\_\_\_\_\_

**Stamp and signature of the authorized electrician:**

# General Terms and Conditions of Service

## I. Scope

Unless otherwise agreed in writing in a specific case, these terms and conditions of service shall apply to service operations, including examining and repairing complaints. All our existing or future legal relationships shall be governed solely by the following terms and conditions of service. Our recognition of any conflicting terms and conditions of the Ordering Party shall be conditional upon our having given our express written consent to their applicability. We hereby expressly object to any terms and conditions of the Ordering Party contained in its General Terms and Conditions of Business or order confirmation. If order confirmations or deliveries are accepted without reservation, this shall not be deemed to constitute recognition of such terms and conditions. Any ancillary agreements or amendments must be confirmed in writing.

## II. Costs

The Ordering Party shall bear the following costs in connection with the service operation:

- De-installation/installation and electrical works (connection / disconnection).
- Transportation, postage and packaging.
- Function testing and troubleshooting including inspection and repair costs.

There shall be no third-party billing.

## III. Obligations / Ordering Party's cooperation

The Ordering Party shall provide free-of-charge assistance to the manufacturer in carrying out the service operation.

In the case of a warranty claim the manufacturer shall provide the required replacement parts to the Ordering Party free of charge.

## IV. Service visit by the manufacturer

In the event that it is essential that a manufacturer employee carry out the service operation on site, this must be agreed in advance. Where the main reason for the service call is not the fault of the manufacturer, any costs incurred shall be recharged to the Ordering Party after the service visit and shall be paid as per agreed payment terms.

## V. Liability

The manufacturer shall assume liability in accordance with the currently applicable statutory regulations. The packaging for all of our products is designed for the shipping of individually packed goods (pallet). We expressly

point out that our packaging is not suitable for individual shipments via parcel post. The manufacturer shall accept no liability for damage incurred as a result of improper packaging in an individual shipment.

## VI. Manufacturer's Guarantee

The manufacturer's guarantee shall apply only in the event that installation, operation and maintenance have been carried out in accordance with the manufacturer's specifications contained in the installation instructions and instructions for use.

- The guarantee period shall commence from the date on which proof of purchase is provided and shall be limited, in principle, to 24 months.
- Guarantee services shall be performed only if the original proof of purchase relating to the equipment can be presented.
- Any and all guarantee claims shall become void if modifications are made to the equipment without the manufacturer's express consent.
- Any guarantee claim shall likewise become void in the case of defects that arise due to repairs or interventions made by unauthorized persons or due to improper use.
- In the case of guarantee claims, the serial and article numbers must be indicated together with the product name and a meaningful description of the fault.
- This guarantee shall cover defective equipment parts, with the exception of usual wear parts. Wear parts are, among others, lamps, glass parts, heating elements and sauna stones.
- Only original replacement parts may be used within the warranty.
- Service visits by outside companies shall require a written order to be issued by our service department.
- The equipment in question shall be sent to our service department by the Ordering Party and at its expense.
- Electrical installation and connection works in the event of service or replacement shall be carried out at the Customer's expense and shall not be borne by the manufacturer.

Complaints in respect of our products shall be reported to the responsible authorized dealer and shall be exclusively handled via the latter.

The manufacturers General Terms and Conditions of Business, which can be found at [www.eos-sauna.com/agb](http://www.eos-sauna.com/agb), shall apply in addition to the foregoing terms and conditions of service.

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